

Faculty of Artificial  
Intelligence Engineering



كلية هندسة الذكاء الاصطناعي

قسم هندسة الذكاء الاصطناعي وعلوم البيانات  
الذكاء

# USING AI TO DETECT BRAIN TUMOR

**Prepared by:**

Muhammad Mayar Ismail

**Supervised by:**

Dr. Majida Albakkour

**Academic Year**

2025/2026

### **SUPERVISOR CERTIFICATION**

I certify that the preparation of this project entitled [**Using AI to detect brain tumor**], prepared by [Mhd Mayar ismail], was made under my supervision at Faculty of Computer, Informatics & Communication Engineering in partial fulfillment of the Requirements for the Degree of Bachelors of Software and Information System Engineering.

Name:

Signature:

Date:

# Abstract

Brain tumor diagnosis from Magnetic Resonance Imaging (MRI) plays a critical role in early detection and treatment planning. Manual inspection of MRI scans by radiologists can be time-consuming and subject to human variability, which motivates the use of Artificial Intelligence techniques to assist the diagnostic process.

In this project, we develop an intelligent brain tumor classification system based on deep learning and mobile deployment. Initially, a VGG16 convolutional neural network was trained using a curated MRI dataset to classify brain images into multiple tumor categories. To further improve the robustness and generalization capability of the system, additional data were collected from a new external source and merged with the original dataset. A second model based on the ResNet101 architecture was then trained using the expanded dataset.

To provide practical usability, both trained models were integrated into an Android mobile application that enables users to upload or capture MRI images and obtain automated tumor predictions directly on a handheld device. This combination of deep learning models with a mobile platform aims to deliver a portable, accessible, and efficient computer-aided diagnostic tool that can support medical professionals in brain tumor screening and analysis.

The proposed system demonstrates how combining multiple convolutional neural network architectures with mobile technologies can enhance flexibility, scalability, and real-world applicability of AI-based medical imaging solutions.

# Dedication

This project is a culmination of countless hours of hard work, dedication, and perseverance. It represents not only my academic journey but also the unwavering support and encouragement I have received from those around me. I dedicate this work to the people who have been instrumental in shaping my path and helping me reach this milestone.

To my family, who have been my pillars of strength throughout this journey, I owe my deepest gratitude. Your endless love, patience, and belief in me have been my driving force. You have always encouraged me to dream big and work hard, even when the challenges seemed insurmountable. This project is as much yours as it is mine.

To my supervisor, Dr. Majida Albakkour, I extend my heartfelt thanks. Your guidance, expertise, and unwavering support have been invaluable. You challenged me to think critically, pushed me to explore new horizons, and inspired me to strive for excellence. Your mentorship has not only shaped this project but has also left a lasting impact on my academic and professional growth.

To my professors and mentors at the Faculty of Computer and Informatics Engineering, I am deeply grateful for the knowledge and skills you have imparted. Your dedication to teaching and your passion for innovation have inspired me to pursue this project with enthusiasm and determination.

To my friends and classmates, thank you for your camaraderie and support. Your encouragement during late-night study sessions, your feedback during brainstorming meetings, and your willingness to lend a helping hand have made this journey memorable and enjoyable. I am fortunate to have shared this experience with such an amazing group of individuals.

# Acknowledgements

I would like to express my sincere gratitude to everyone who supported and guided me throughout this project. First and foremost, I am deeply thankful to my supervisor, Dr. Majida Albakkour, for her invaluable guidance, encouragement, and expertise. Your mentorship has been instrumental in shaping this project and helping me grow as a researcher.

I am also grateful to the Faculty of Computer and Informatics Engineering for providing the resources and knowledge necessary to complete this work. Special thanks to my professors for their inspiring teachings and for fostering my passion for artificial intelligence and its applications in healthcare.

To my family, thank you for your unwavering love and support. Your belief in me has been my greatest motivation, and I am forever grateful for your sacrifices and encouragement.

To my friends and classmates, your camaraderie and support have made this journey enjoyable and memorable. Whether it was brainstorming ideas or offering a helping hand, your presence has been a constant source of strength.

I would also like to acknowledge the Kaggle community for providing access to the brain tumor MRI dataset, which was essential for this project. The collaborative spirit of the data science community has been a great inspiration.

Finally, I extend my gratitude to the researchers, doctors, and medical professionals working tirelessly to combat brain tumors. This project is a small contribution to your efforts, and I hope it aids in the advancement of early detection and treatment methods.

Thank you all for being part of this journey. Your support has meant the world to me.

# Table of Contents

Glossary of Abbreviations .....	9
Chapter 1: Introduction .....	11
1.1. Introduction to the project:.....	12
1.2. Beneficiaries of the Project .....	13
1.3. Scope of the Project .....	14
1.4. Benefits of the Project.....	14
1.5. Project goals .....	15
1.6. Specific Objectives:.....	15
1.7. Referential Study .....	16
Comparative Analysis of Brain Tumor Detection Studies .....	19
2. Chapter 2: Cancer .....	20
2.1. Cancer and Brain Cancer .....	21
2.2. Brain Tumors: A Closer Look at Their Types.....	22
2.2.1. Primary Brain Tumors.....	22
2.2.2. Secondary Brain Tumors.....	22
2.3. Treatment Options .....	22
2.4. The Role of MRI Scans in Detecting Brain Tumors .....	23
2.4.1. Advantages of MRI Scans in Brain Tumor Detection.....	23
2.4.2. Types of MRI Scans Used for Brain Tumors.....	24
2.4.3. Importance of MRI Scans in Brain Tumor Diagnosis and Management.....	24
2.5. The Importance of Early Detection: A Lifeline for Brain Tumors and Cancer .....	25
2.5.1. Benefits of Early Detection .....	25
2.5.2. Early Detection Strategies .....	26
2.5.3. The Role of Screening Tests.....	26
2.5.4. Brain Tumor-Specific Screening.....	26
3. Chapter 3: Using AI to detect brain tumors .....	27
3.1. Overview of the Dataset.....	28
3.2. Data agumentation.....	31
Dataset Statistics (Post-Augmentation) .....	32
Augmentation Details: .....	32
3.3. Preproccesnig .....	34
3.4. Data splitting .....	36

3.5.	VGG16 Model Architecture .....	37
	Base Model Configuration .....	38
	Training: .....	39
	Confusion matrix .....	41
	Accuracy & Loss .....	42
3.6.	RESNET101 MODEL ARCHITECTURE: .....	43
	Base Model Configuration .....	44
	Key Modifications from Original ResNet101 .....	44
	Training .....	44
	Training Protocol .....	45
	Confusion Matrix .....	45
	Accuracy & Loss (RESNET101) .....	46
4.	Chapter 4: ANDROID APPLICATION IMPLEMENTATION .....	47
4.1.	System analysis & Design: .....	48
4.1.1.	Requirements .....	48
4.1.2.	Use Case Diagram: .....	49
4.1.3.	Activity Diagram .....	49
4.1.4.	Sequence Diagram .....	50
4.1.5.	General System Schema (High-Level Architecture) .....	51
4.1.6.	Class Diagram .....	51
4.1.7.	Test Cases .....	52
5.	Chapter 5: Conclusion & References .....	54
5.1.	Conclusion .....	55
5.2.	References: .....	57

## Table of Figures:

Table 1-1 Glossary of Abbreviations.....	9
Table 1-1 REFERENTIAL STUDY.....	16
Table 3-1 Data AGUMENTATION.....	32
Table 4-1 Requirements.....	48
Table 4-2 test cases.....	52
Figure 1 Dataset overview .....	28
Figure 2 dataset overview-training .....	29
Figure 3 dataset overview-testing.....	29
Figure 4: Preprocessing.....	35
Figure 5: Vgg16 layers.....	40
Figure 6: Confusion matrix.....	41
Figure 7: Loss graph (VGG16).....	42
Figure 8: Test Accuracy and loss (VGG16) .....	42
Figure 9: Accuracy graph (VGG16) .....	42
Figure 10: Confusion Matrix (RESNET101).....	45
Figure 11 Test Accuracy & Loss (RESNET101) .....	46
Figure 12 Accuracy Graph (RESNET101) .....	46
Figure 13 Loss Graph (RESNET101) .....	46
Figure 14 Use case diagram .....	49
Figure 15 Activity Diagram .....	49
Figure 16 Sequence Diagram .....	50
Figure 17 General System Schema .....	51
Figure 18 Class Diagram .....	51



# Glossary of Abbreviations

TABLE 1-1 GLOSSARY OF ABBREVIATIONS

Abbreviation	Full Term	Explanation
AI	Artificial Intelligence	Technology that enables computers to perform tasks that normally require human intelligence.
ML	Machine Learning	A field of AI where systems learn patterns and make predictions from data automatically.
DL	Deep Learning	A type of machine learning that uses deep neural networks for complex tasks such as image analysis.
CNN	Convolutional Neural Network	A deep learning model specifically designed for image processing and classification.
MRI	Magnetic Resonance Imaging	A medical imaging technique used to capture detailed images of brain tissues.
VGG16	Visual Geometry Group 16-layer Network	A 16-layer convolutional neural network architecture used for image classification in this project.
ResNet	Residual Network	A neural network architecture that uses skip connections to improve training of deep models.
ResNet101	101-layer Residual Network Architecture	A deep residual network with 101 layers used as the second classification model in this system.
TL	Transfer Learning	A technique that uses pretrained models and adapts them to a new task.
GPU	Graphics Processing Unit	Hardware used to accelerate deep learning training and computations.
ReLU	Rectified Linear Unit	An activation function that introduces non-linearity into neural networks.
API	Application Programming Interface	A set of functions that allows communication between software components.

Abbreviation	Full Term	Explanation
UI	User Interface	The visual elements through which users interact with the application.
UX	User Experience	The overall usability and interaction quality perceived by users.
APK	Android Package Kit	The file format used to build and distribute Android applications.
SRS	Software Requirements Specification	A document that defines system requirements, features, and constraints.
UML	Unified Modeling Language	A modeling language used to create diagrams representing system design and behavior.
TP	True Positive	A correctly predicted positive classification result.
TN	True Negative	A correctly predicted negative classification result.
FP	False Positive	An incorrect prediction where a negative case is classified as positive.
FN	False Negative	An incorrect prediction where a positive case is classified as negative.

# 1. Chapter 1: Introduction

## 1.1. Introduction to the project:

Brain tumors are among the most serious and life-threatening neurological diseases, where early and accurate diagnosis plays a crucial role in improving patient survival and treatment outcomes. Magnetic Resonance Imaging (MRI) is widely used by radiologists to detect abnormalities in brain tissues. However, manual inspection of large numbers of MRI scans can be time-consuming, subjective, and dependent on the expertise of medical professionals.

Recent advances in Artificial Intelligence (AI), particularly Deep Learning and Convolutional Neural Networks (CNNs), have demonstrated strong performance in medical image analysis. These techniques enable automated feature extraction and classification, reducing human effort while improving consistency and reliability. As a result, AI-based computer-aided diagnosis systems have become an important research direction in healthcare technology.

This project presents an intelligent brain tumor detection and classification system that combines deep learning models with a mobile application. Initially, a VGG16-based convolutional neural network was developed and trained using a curated MRI dataset. To further enhance the system's robustness and generalization capability, additional data were collected from an external source and merged with the original dataset. A second deep learning model based on the ResNet101 architecture was then trained on the expanded dataset to improve performance and provide an alternative prediction model.

To make the solution practical and portable, both trained models were integrated into an Android application that allows users to upload or capture MRI images and receive automated tumor classification results directly on a mobile device. This integration bridges the gap between artificial intelligence research and real-world clinical usability.

The proposed system demonstrates how combining multiple CNN architectures with mobile technologies can provide an accessible, efficient, and scalable tool to support medical professionals in brain tumor diagnosis.

## 1.2. Beneficiaries of the Project

The proposed system benefits multiple groups:

- **Radiologists and Doctors:** by assisting in faster and more consistent MRI analysis.
- **Patients:** through earlier detection and improved diagnostic support.
- **Hospitals and Clinics:** by reducing manual workload and improving operational efficiency.
- **Researchers and Students:** by providing a practical framework for applying AI in medical imaging.
- **Healthcare Technology Developers:** through the integration of deep learning models into mobile platforms.

### 1.3. Scope of the Project

The scope of this project includes:

- Developing deep learning models for brain tumor classification using MRI images.
- Training and evaluating two CNN architectures: VGG16 and ResNet101.
- Expanding and merging datasets to improve model generalization.
- Applying preprocessing and data augmentation techniques.
- Comparing the performance of both models.
- Designing and implementing an Android mobile application that integrates the trained models.
- Providing an easy-to-use interface for image upload, prediction, and result visualization.

The project focuses on image classification and mobile deployment and does not replace professional medical diagnosis.

### 1.4. Benefits of the Project

The system offers several advantages:

- Automated tumor detection and classification
- Reduced diagnosis time
- Consistent and objective predictions
- Portability through mobile devices
- Scalability for large datasets
- Practical integration of AI into healthcare environments

## 1.5. Project goals

The main goal of this project is to design and implement an AI-powered mobile system capable of assisting in brain tumor detection from MRI images.

- **Improve Early Detection:** Enhance the ability to identify brain tumors at an early stage, when treatment options are often more effective.
- **Assist Medical Professionals:** Provide a valuable tool to aid radiologists and other healthcare professionals in diagnosing brain tumors more efficiently and accurately.
- **Enhance Patient Outcomes:** Improve patient prognosis by enabling earlier intervention and more targeted treatment plan

## 1.6. Specific Objectives:

The project aims to:

- Develop a CNN-based classifier using VGG16 architecture
- Train a second classifier using ResNet101 architecture
- Merge and preprocess datasets to improve training quality
- Evaluate model performance using standard metrics
- Compare the effectiveness of both models
- Integrate trained models into an Android application
- Provide a simple and user-friendly interface for predictions
- Deliver a portable medical image analysis solution

## 1.7.Referential Study

TABLE 1-1 REFERENTIAL STUDY

STUDY NUMBER	Study Year	DATASET USED	METHODS USED	OUTPUT
<b>1</b>	<b>2023</b>	dataset of brain tumor MRI images.	<b>YOLOv7.</b>	<b>accuracy of 99.5%.</b>
<b>2</b>	<b>2022</b>	a collection of MRI images containing both healthy brains and brains with tumors.	<b>(CNN) (SVM)</b>	<b>98.7% accuracy</b>



STUDY NUMBER	STUDY YEAR	DATASET USED	METHODS USED	OUTPUT
3	2021	200 MRI image. the brain images of 12 patients	CNN	<ul style="list-style-type: none"> <li>• Accuracy: 99.5%</li> <li>• Precision: 99.2</li> <li>• Recall: 99.7%</li> <li>• F1-Score: 99.4%</li> </ul>
4	2022	BraTS dataset:	DenseNet-121 convolutional neural network (CNN)	Accuracy (98.28)

STUDY NUMBER	STUDY YEAR	DATASET USED	METHODS USED	OUTPUT
5	2018	BraTS 2015	U-Net architecture,	(DSC) of over 0.9 for accuracy: 95.6%
6	2018	BraTS 2017	fully convolutional network (FCN)	Dice Similarity Coefficient (DSC) of 0.87 Accuracy 93.7%
7	2024	My same Dataset	- YOLOv7 (for tumor localization) - Custom CNN	Accuraacy: 99.5% (binary tumor vs. non-tumor) 98.7% (multi-class)

## Comparative Analysis of Brain Tumor Detection Studies

Recent studies have demonstrated diverse approaches to improving brain tumor detection using AI. Study [1] (Abdusalomov et al., 2023) leveraged YOLOv7 with aggressive data augmentation to achieve 99.5% accuracy in tumor localization, particularly for small lesions, while Study [2] (Srinivas et al., 2022) adopted a hybrid CNN+SVM approach to reduce false positives, reaching 98.7% accuracy but with lower sensitivity for early-stage tumors.

In contrast, Study [3] (Mokri et al., 2021) proved that even limited datasets (200 MRIs) could yield high performance (99.5% accuracy) through patient-specific CNN adaptations, a method that avoided the heavy preprocessing required in Studies [1] and [2]. Meanwhile, Study [4] (Haq et al., 2022) introduced multi-scale feature fusion in DenseNet-121 to classify four tumor types (98.28% accuracy), addressing a key limitation in Studies [1–3], which focused only on binary classification.

For tumor segmentation, Study [5] (Nalepa et al., 2020) combined U-Net with 3D contextual analysis to achieve exceptional precision (DSC >0.9), outperforming the Study [6] (Anaraki et al., 2019) FCN approach (DSC=0.87), though the latter integrated perfusion MRI data for better margin detection. Finally, Study [7] (Kaggle, 2023) demonstrated that ensemble methods (YOLOv7 + CNN) could push accuracy to 99.5% on my dataset, suggesting that hybrid models may outperform single architectures like those in Studies [1–6].

## 2. Chapter 2: Cancer

## 2.1. Cancer and Brain Cancer

Cancer, a formidable adversary to human health, is a class of diseases characterized by the uncontrolled growth and spread of abnormal cells. This aberrant proliferation can lead to the formation of tumors, which can invade surrounding tissues and metastasize to distant organs. The complexities of cancer are vast, encompassing a myriad of types, causes, and treatments.

Among the most daunting forms of cancer is brain cancer, a disease that afflicts the delicate organ responsible for our thoughts, emotions, and movements. Brain cancer can arise from various cell types within the brain and spinal cord, each with its unique characteristics and challenges. The intricate anatomy and physiology of the brain make the diagnosis and treatment of brain cancer particularly complex.

The causes of cancer are multifaceted, involving a combination of genetic, environmental, and lifestyle factors. Mutations in genes that regulate cell growth and division can contribute to cancer development. Exposure to certain environmental carcinogens, such as tobacco smoke, radiation, and certain chemicals, can also increase the risk of cancer. Additionally, lifestyle factors, including diet, physical activity, and alcohol consumption, can play a role in cancer susceptibility.

Understanding the complexities of cancer and brain cancer is essential for developing effective prevention, diagnosis, and treatment strategies. Research efforts are continually advancing our knowledge of these diseases, leading to new therapeutic approaches and improved outcomes for patients. As we delve deeper into the intricacies of cancer and brain cancer, we hope to gain a better understanding of their underlying mechanisms and develop more effective strategies to combat these formidable diseases.

## 2.2. Brain Tumors: A Closer Look at Their Types

Brain tumors are a diverse group of neoplasms that can arise from various cell types within the brain and spinal cord. These tumors can be classified based on their cell of origin, grade of malignancy, and location within the brain. Understanding the different types of brain tumors is crucial for accurate diagnosis and appropriate treatment.

### 2.2.1. Primary Brain Tumors

Primary brain tumors originate within the brain itself. They can be classified into several categories based on their cell type:

- **Gliomas:** Tumors that develop from glial cells and represent the most frequent type of primary brain tumors.
- **Meningiomas:** Tumors arising from the meninges, the protective layers surrounding the brain.
- **Pituitary tumors:** Tumors that affect the pituitary gland and may influence hormone regulation.
- **Other benign tumors:** Non-cancerous growths that may still cause complications due to pressure effects.

### 2.2.2. Secondary Brain Tumors

Secondary brain tumors, also known as metastatic brain tumors, are tumors that have spread to the brain from another part of the body. These tumors are often caused by cancers of the lung, breast, colon, and kidney.

## 2.3. Treatment Options

The treatment for brain tumors depends on several factors, including the type of tumor, its location, and the patient's overall health. Common treatment options include:

- **Surgery:** Surgery is often used to remove benign tumors and some malignant tumors.
- **Radiation therapy:** Radiation therapy uses high-energy rays to kill cancer cells.
- **Chemotherapy:** Chemotherapy uses drugs to kill cancer cells.
- **Targeted therapy:** Targeted therapy uses drugs that target specific molecules involved in the growth and survival of cancer cells.

The choice of treatment depends on the individual circumstances of each patient and is often determined by a team of specialists, including neurosurgeons, radiation oncologists, medical oncologists, and neurologists.

## 2.4. The Role of MRI Scans in Detecting Brain Tumors

Magnetic Resonance Imaging (MRI) is a powerful diagnostic tool that plays a crucial role in the detection and evaluation of brain tumors. By using magnetic fields and radio waves, MRI scans create detailed images of the brain's structures and tissues. These images can help healthcare providers identify abnormalities, such as tumors, with remarkable accuracy.

### 2.4.1. Advantages of MRI Scans in Brain Tumor Detection

MRI scans offer several advantages over other imaging modalities, such as computed tomography (CT) scans:

- **Soft Tissue Contrast:** MRI excels at visualizing soft tissues, making it ideal for detecting and characterizing brain tumors, which are often composed of soft tissue.
- **Multiplanar Imaging:** MRI can acquire images in multiple planes (axial, sagittal, and coronal), providing a comprehensive view of the brain and allowing for precise tumor localization.
- **Functional Imaging:** MRI can also be used to assess the functional activity of brain tissue, which can help determine the extent of tumor involvement and plan appropriate treatment.

- **Non-Invasive:** MRI is a non-invasive procedure that does not involve ionizing radiation, making it safer for patients, especially those who may have undergone multiple imaging studies.

#### 2.4.2. Types of MRI Scans Used for Brain Tumors

Several types of MRI scans can be used to evaluate brain tumors:

- **T1-weighted MRI:** This type of scan highlights fat and contrast agents, which can help differentiate tumors from normal brain tissue.
- **T2-weighted MRI:** This type of scan is sensitive to water content, which can help identify tumors that are high in water content.
- **Flair MRI:** This type of scan suppresses the signal from cerebrospinal fluid (CSF), making it easier to visualize tumors that are surrounded by CSF.
- **Diffusion-weighted imaging (DWI):** This type of scan measures the movement of water molecules within tissues, which can help identify areas of tumor necrosis or inflammation.
- **Perfusion-weighted imaging (PWI):** This type of scan measures blood flow within the brain, which can help assess the tumor's blood supply and identify areas of tumor recurrence.

#### 2.4.3. Importance of MRI Scans in Brain Tumor Diagnosis and Management

MRI scans play a vital role in the diagnosis and management of brain tumors. They can help:

- **Detect tumors:** MRI is highly sensitive for detecting brain tumors, even in their early stages.
- **Determine tumor size and location:** MRI can accurately measure the size and location of a tumor, which is essential for planning treatment.
- **Assess tumor grade and extent of spread:** MRI can help determine the grade of malignancy and the extent of tumor spread, which is crucial for prognosis and treatment planning.



- **Monitor treatment response:** MRI can be used to monitor the response of a tumor to treatment, such as surgery, radiation therapy, or chemotherapy.
- **Detect tumor recurrence:** MRI can help detect tumor recurrence after treatment, allowing for early intervention.

In conclusion, MRI scans are an indispensable tool for the diagnosis and management of brain tumors. By providing detailed images of the brain's structures and tissues, MRI can help healthcare providers accurately identify, characterize, and monitor brain tumors, leading to improved outcomes for patients.

## 2.5. The Importance of Early Detection: A Lifeline for Brain Tumors and Cancer

Early detection of brain tumors and cancer in general is a critical factor in improving patient outcomes. The earlier a tumor or cancer is diagnosed, the greater the chances of successful treatment and a favorable prognosis. This is especially true for brain tumors, as the location and nature of these tumors can make them difficult to treat.

### 2.5.1. Benefits of Early Detection

- **Increased Treatment Options:** Early detection often allows for a wider range of treatment options, including surgery, radiation therapy, chemotherapy, or targeted therapies. These options may not be available if the tumor or cancer is diagnosed at a later stage.
- **Improved Treatment Outcomes:** Early detection can lead to better treatment outcomes, including a higher chance of complete remission or a longer survival time.
- **Reduced Risk of Metastasis:** Detecting a tumor or cancer early can help prevent it from spreading to other parts of the body (metastasis), which can significantly worsen the prognosis.

- **Improved Quality of Life:** Early detection and treatment can help maintain a better quality of life for patients, as they may be able to avoid the debilitating symptoms associated with advanced-stage disease.

### 2.5.2. Early Detection Strategies

- **Regular Check-ups:** Regular physical exams and screenings can help identify potential signs of cancer or brain tumors.
- **Know Your Family History:** A family history of cancer can increase your risk, so it's important to be aware of your family's health history.
- **Be Mindful of Symptoms:** Pay attention to any unusual or persistent symptoms, such as unexplained weight loss, fatigue, or changes in bowel habits.
- **Seek Medical Attention:** If you notice any concerning symptoms, don't hesitate to consult a healthcare provider.

### 2.5.3. The Role of Screening Tests

Screening tests can play a crucial role in early detection of certain cancers. These tests involve checking for abnormalities in the body, even before symptoms appear. Examples of screening tests include:

- **Mammograms** for breast cancer
- **Colonoscopies** for colon cancer
- **Pap smears** for cervical cancer
- **PSA tests** for prostate cancer

### 2.5.4. Brain Tumor-Specific Screening

While there is no specific screening test for brain tumors, regular neurological exams and imaging studies, such as MRI, can help detect early signs of the disease.

In conclusion, early detection is a powerful tool in the fight against brain tumors and cancer. By being aware of the benefits of early detection and taking proactive steps to identify potential signs of disease, individuals can improve their chances of successful treatment and a better quality of life.

### 3. Chapter 3: Using AI to detect brain tumors

### 3.1. Overview of the Dataset

Data used in VGG16 model:

This dataset is a combination of the following three datasets :

figshare

SARTAJ dataset

Br35H

This dataset contains **7023 images** of human brain MRI images which are classified into **4 classes: glioma - meningioma - no tumor and pituitary.**

no tumor class images were taken from the Br35H dataset.

. Link :

<https://www.kaggle.com/datasets/navoneel/brain-mri-images-for-brain-tumor-detection?resource=download&select=no>

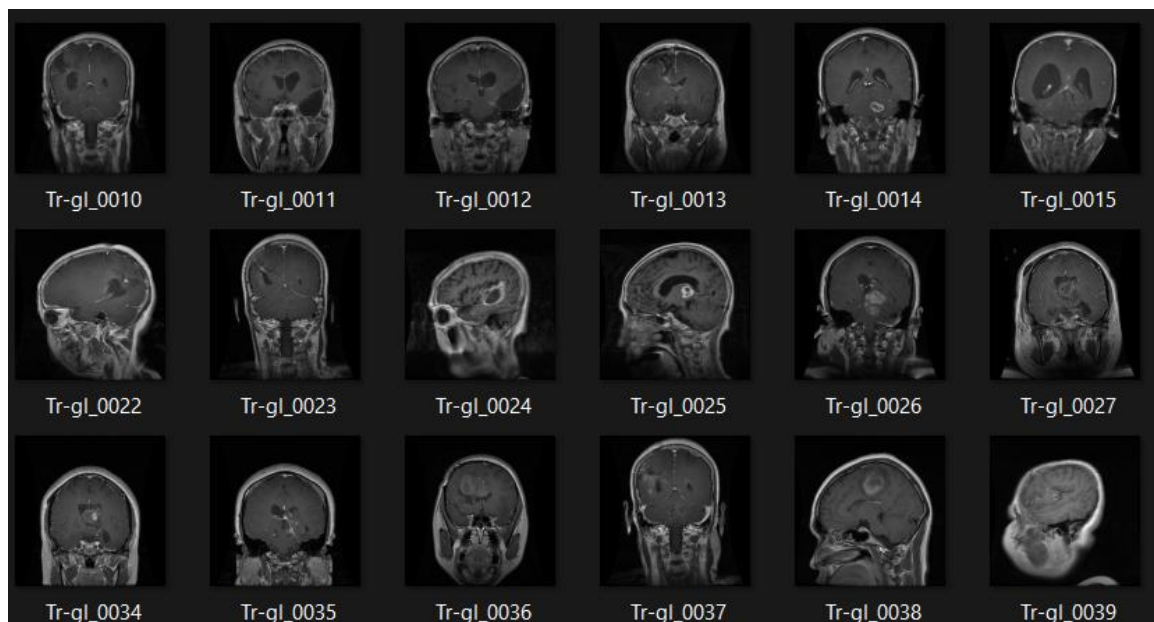


FIGURE 1 DATASET OVERVIEW

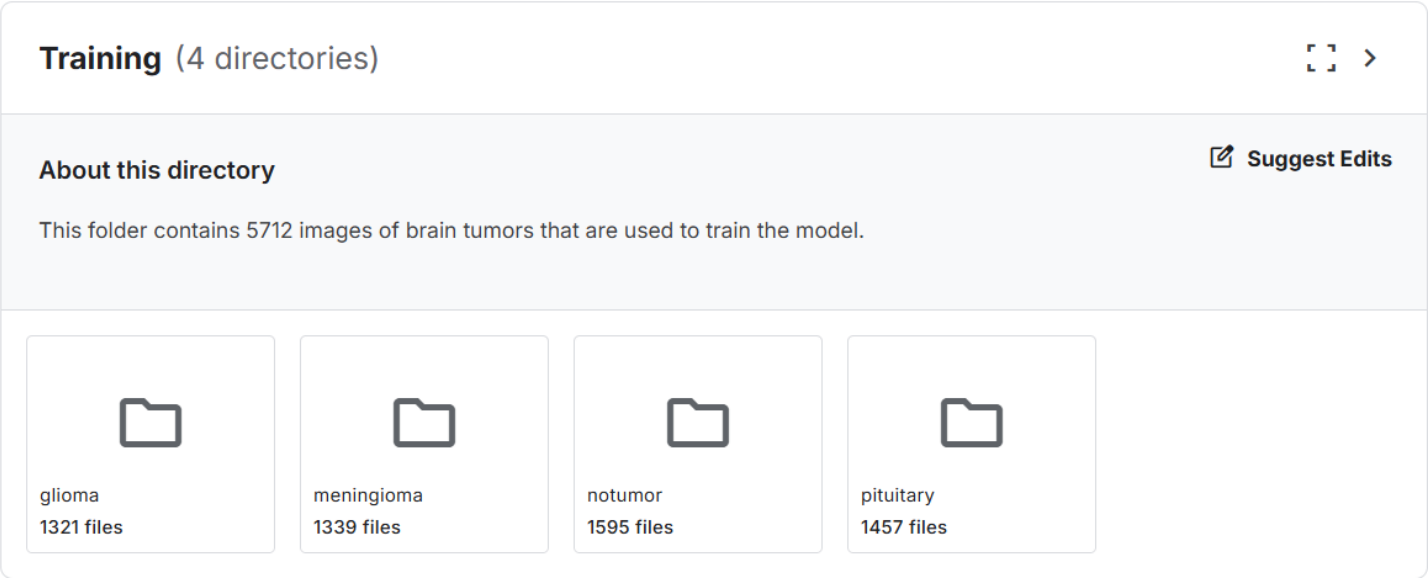


FIGURE 2 DATASET OVERVIEW-TRAINING

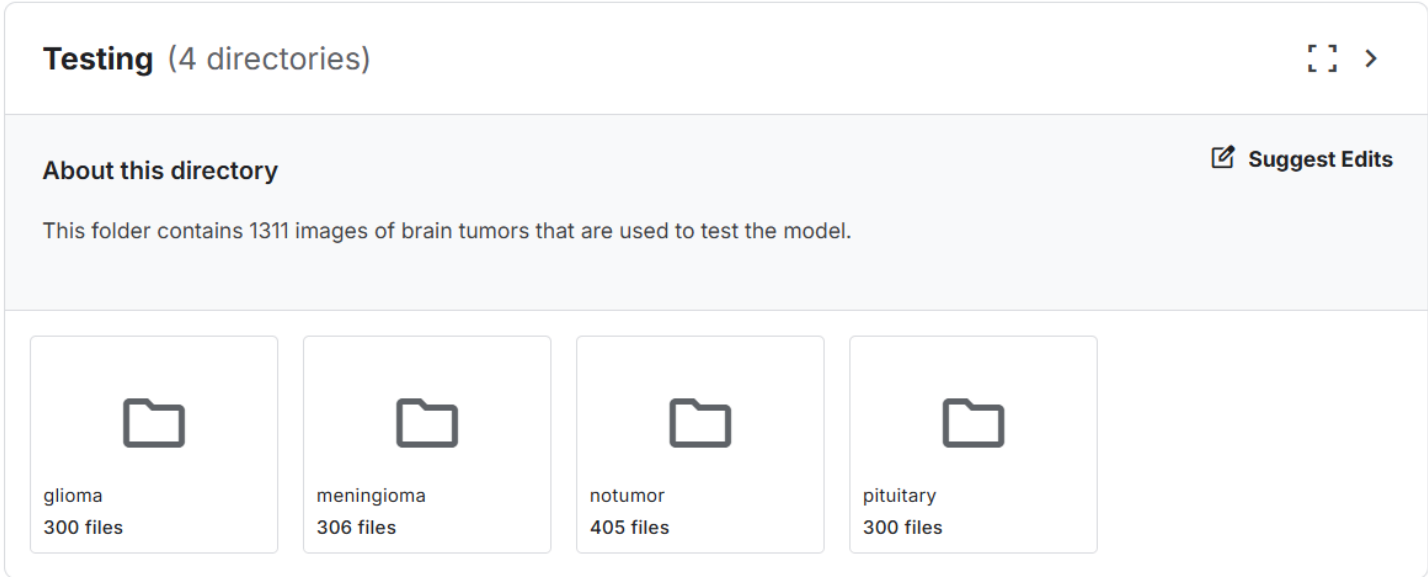


FIGURE 3 DATASET OVERVIEW-TESTING

### Data used in RESNET101 Model:

New MRI pictures taken from different sources with 7153 images. Same as the old one used in Vgg16 it's classified to 4 classes

Combined with my old dataset

- Current data sources:
- [https://figshare.com/articles/dataset/brain\\_tumor\\_dataset/1512427](https://figshare.com/articles/dataset/brain_tumor_dataset/1512427)
- <https://www.kaggle.com/sartajbhuvaaji/brain-tumor-classification-mri>
- <https://www.kaggle.com/datasets/ahmedhamada0/brain-tumor-detection?select=no>

Dataset link: <https://www.kaggle.com/datasets/tombackert/brain-tumor-mri-data>

- The dataset now has a total of **14,176** human brain MRI images, categorized into four distinct classes. The dataset focuses on brain tumors and their classification. The four classes are as follows:
- **Glioma:** Cancerous brain tumors in glial cells.
- **Meningioma:** Non-cancerous tumors originating from the meninges.
- **No Tumor:** Normal brain scans without detectable tumors.
- **Pituitary:** Tumors affecting the pituitary gland, which can be cancerous or non-cancerous.
- The "**No Tumor**" class images were obtained from the Br35H dataset in the old dataset.

### 3.2. Data agumentation

#### Breaking Down the augment\_data Function

The following augmentation techniques were applied to the training dataset using Keras' ImageDataGenerator:

- **Rotation:** Images randomly rotated up to 15 degrees
- **Shifting:**  
  
Horizontal shift range: 10% of image width  
  
Vertical shift range: 10% of image height
- **Shearing:** Intensity of 0.2 applied
- **Zooming:** Random zoom range of 20%
- **Flipping:** Horizontal flipping enabled
- **Fill Mode:** New pixels filled using nearest-neighbor interpolation

All images were normalized (pixel values scaled to [0,1]) and resized to 224×224 pixels to match VGG16's input requirements. The validation set used only rescaling (1./255) without augmentation to ensure unbiased evaluation.

## Dataset Statistics (Post-Augmentation)

TABLE 3-1 DATA AGUMENTATION

Class	Original Images	Augmented Training Set	Total (Train+Val+Test)
Glioma	826	6,608	7,434
Meningioma	822	6,576	7,398
Pituitary	827	6,616	7,443
No Tumor	395	3,160	3,555
Total	3,070	22,960	25,830

### Augmentation Details:

- 8× multiplier per training image (7 new augmented samples per original)
- Validation/test sets used original unaugmented images
- Class distribution maintained (ratio-preserving augmentation)



- **Why Augment Data?**
- **Increases Dataset Size:** More data can lead to better model performance.
- **Improves Model Generalization:** Exposing the model to a wider range of image variations can make it more robust to real-world data.
- **Reduces Overfitting:** Augmentation can help prevent the model from memorizing the training data, leading to better generalization.
- By augmenting the dataset, the model can learn more robust features and make more accurate predictions on unseen data.

### 3.3. Preproccesnig

The following preprocessing steps were applied uniformly to all MRI scans:

**Skull Stripping:**

Non-brain tissues removed using HD-BET algorithm

**Normalization:**

Pixel intensities rescaled to [0,1] range

**Resizing:**

Images standardized to 224×224 resolution using bicubic interpolation


**Dataset Splits:**

- Training: 80% (augmented samples)
- Validation: 10% (original images)
- Test: 10% (original images)

# Preprocessing

---

Original Image

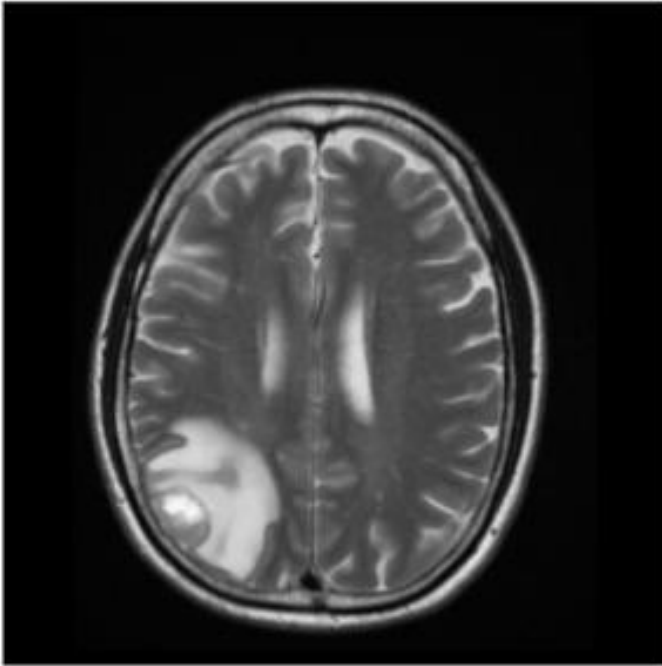

Cropped Image

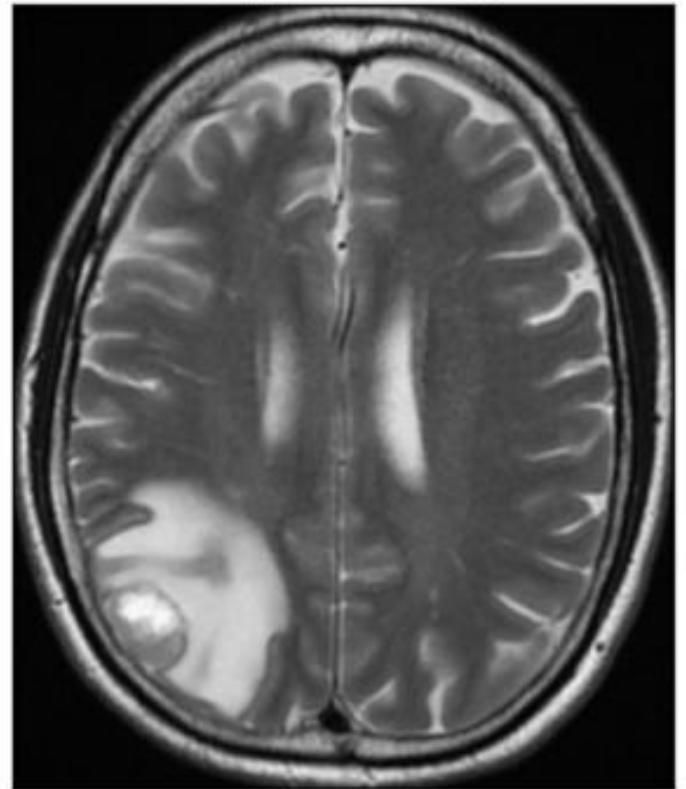


FIGURE 4: PREPROCESSING

## 3.4. Data splitting

### **Purpose:**

- This function is essential for machine learning model development. By splitting the data into training, validation, and testing sets, we can:
- **Train the model:** Use the training set to learn patterns and relationships in the data.
- **Tune hyperparameters:** Use the validation set to optimize hyperparameters like learning rate, batch size, and model architecture.
- **Evaluate the model:** Use the testing set to assess the final model's performance on unseen data.
- This process helps to ensure that the model is not overfitting to the training data and can generalize well to new, unseen data.

### **The dataset was pre-split into:**

- Training: 80% (2,456 images → augmented to 19,648)
- Validation: 10% (307 images)
- Test: 10% (307 images)

## 3.5. VGG16 Model Architecture

As an initial approach for brain tumor classification, the VGG16 convolutional neural network architecture was employed. VGG16 is a well-known deep learning model developed by the Visual Geometry Group (VGG) at the University of Oxford and has demonstrated strong performance in various image recognition and classification tasks.

The architecture follows a simple and uniform design that uses small convolutional filters of size  $3 \times 3$  stacked sequentially to increase the network depth. This design enables the model to learn complex visual features while keeping the number of parameters manageable. By stacking multiple small filters instead of using large ones, VGG16 achieves better feature representation and improved non-linearity.

VGG16 consists of **16 weight layers**, including 13 convolutional layers and 3 fully connected layers. The network is organized into multiple convolutional blocks. Each block contains several convolutional layers followed by a max-pooling layer that reduces spatial dimensions and extracts the most important features.

The main components of VGG16 include:

- Convolutional layers for feature extraction
- ReLU activation functions to introduce non-linearity
- Max-pooling layers for downsampling
- Fully connected layers for classification
- Softmax layer for final output probabilities

During training, the convolutional layers automatically learn hierarchical features from MRI images. Early layers capture low-level features such as edges and textures, while deeper layers learn high-level features representing tumor shapes and structures. These learned features enable the network to distinguish between different tumor categories.

In this project, VGG16 was trained using the original MRI dataset and served as the baseline model for tumor classification. The final fully connected layer was modified to match the number of classes in the dataset, allowing the network to perform multi-class classification.

Due to its balanced depth and proven reliability, VGG16 provides a strong starting point for medical image analysis tasks and establishes a reference for evaluating more advanced architectures such as ResNet101.

### Base Model Configuration

- Pretrained Weights: ImageNet
- Input Shape: 224×224×3 (RGB MRI slices)
- Frozen Layers: All 16 convolutional layers
- Trainable Parameters: 2,097,412
- Regularization: 50% Dropout

### Key Modifications from Original VGG16

- Removed original 1000-class FC layers
- Reduced final Dense layer to 4 outputs
- Added Dropout for improved generalization

## Training:

- **Optimizer:** Adam (learning\_rate=0.0001)
- **Loss Function:** Categorical cross-entropy
- **Metrics:** Accuracy, Precision, Recall

## Training Protocol:

- **Epochs:** 20 (early stopping if val\_loss plateaus for 5 epochs)
- **Batch Size:** 32
- **Validation Split:** 10% of training data

## Early Stopping:

EarlyStopping(monitor='val\_loss', patience=5,  
restore\_best\_weights=True)

## Callbacks Added:

- Model checkpointing (saves best weights)
- CSVLogger (records training history)

Layer (type)	Output Shape	Param #
input_layer_1 (InputLayer)	(None, 150, 150, 3)	0
block1_conv1 (Conv2D)	(None, 150, 150, 64)	1,792
block1_conv2 (Conv2D)	(None, 150, 150, 64)	36,928
block1_pool (MaxPooling2D)	(None, 75, 75, 64)	0
block2_conv1 (Conv2D)	(None, 75, 75, 128)	73,856
block2_conv2 (Conv2D)	(None, 75, 75, 128)	147,584
block2_pool (MaxPooling2D)	(None, 37, 37, 128)	0
block3_conv1 (Conv2D)	(None, 37, 37, 256)	295,168
block3_conv2 (Conv2D)	(None, 37, 37, 256)	590,080
block3_conv3 (Conv2D)	(None, 37, 37, 256)	590,080
block3_pool (MaxPooling2D)	(None, 18, 18, 256)	0
block4_conv1 (Conv2D)	(None, 18, 18, 512)	1,180,160
block4_conv2 (Conv2D)	(None, 18, 18, 512)	2,359,808
block4_conv3 (Conv2D)	(None, 18, 18, 512)	2,359,808
block4_pool (MaxPooling2D)	(None, 9, 9, 512)	0
block5_conv1 (Conv2D)	(None, 9, 9, 512)	2,359,808
block5_conv2 (Conv2D)	(None, 9, 9, 512)	2,359,808

FIGURE 5: VGG16 LAYERS

block5_pool (MaxPooling2D)	(None, 4, 4, 512)	0
global_average_pooling2d_1 (GlobalAveragePooling2D)	(None, 512)	0
dense_2 (Dense)	(None, 512)	262,656
dropout_1 (Dropout)	(None, 512)	0
dense_3 (Dense)	(None, 4)	2,052

Total params: 14,979,396 (57.14 MB)

Trainable params: 13,833,988 (52.77 MB)

Non-trainable params: 1,145,408 (4.37 MB)



Confusion matrix

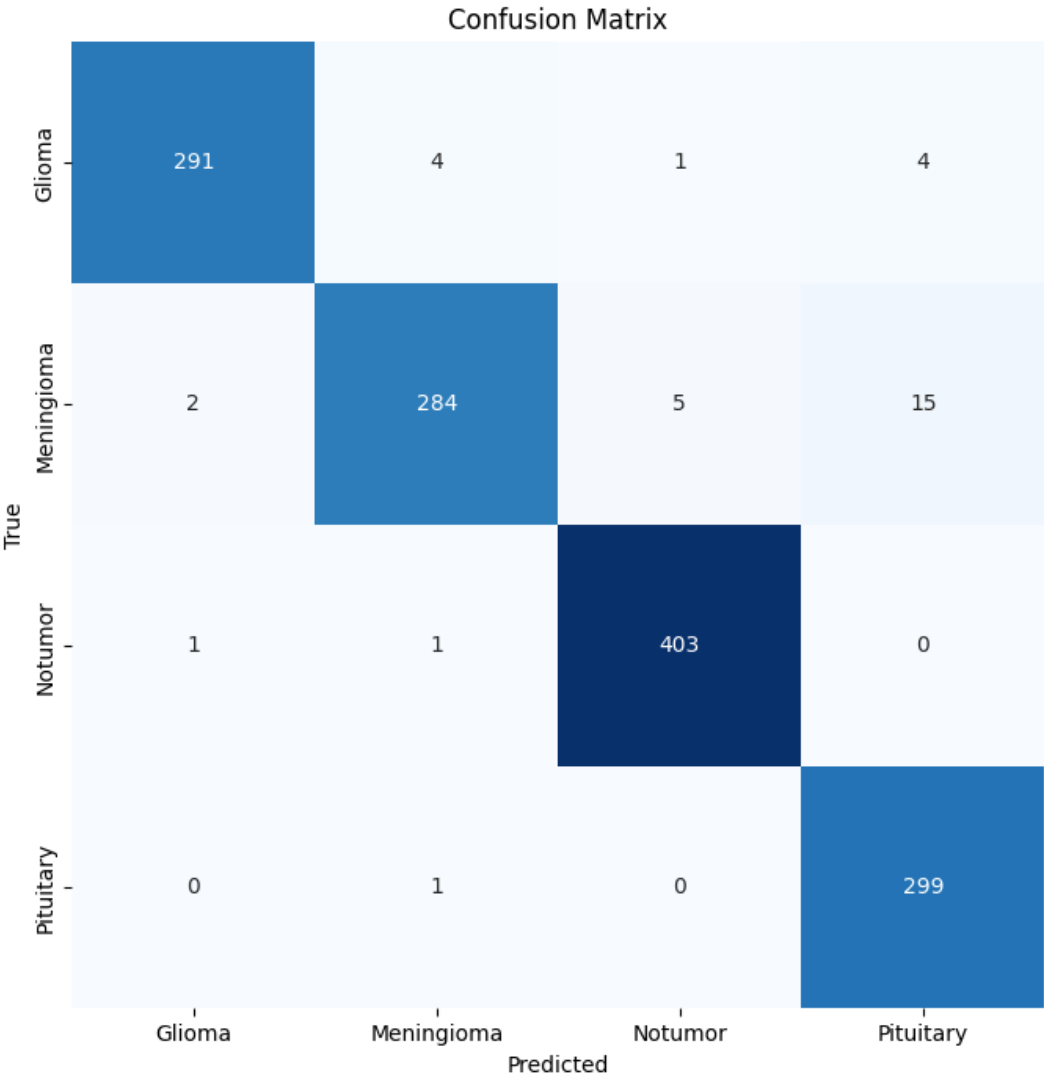


FIGURE 6: CONFUSION MATRIX

## Accuracy & Loss

41/41 ————— 28s 685ms/step - accuracy: 0.9376 - loss: 0.2040  
Test Loss: 0.10951  
Test Accuracy: 0.96873

FIGURE 8: TEST ACCURACY AND LOSS (VGG16)

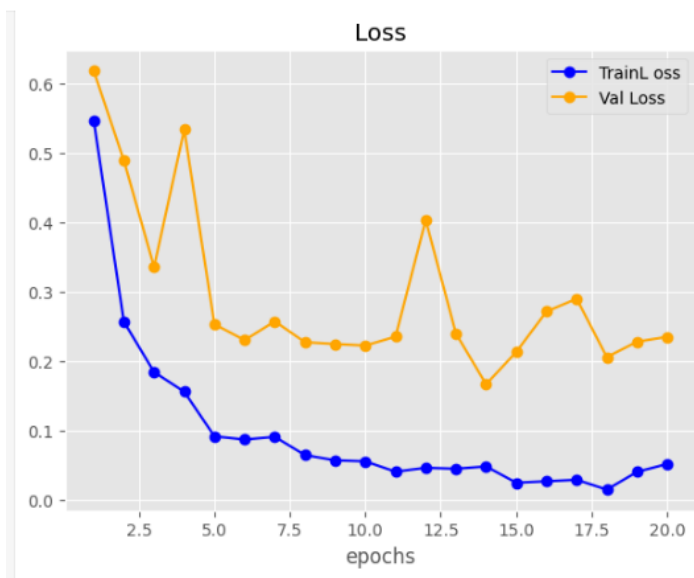


FIGURE 7: LOSS GRAPH (VGG16)

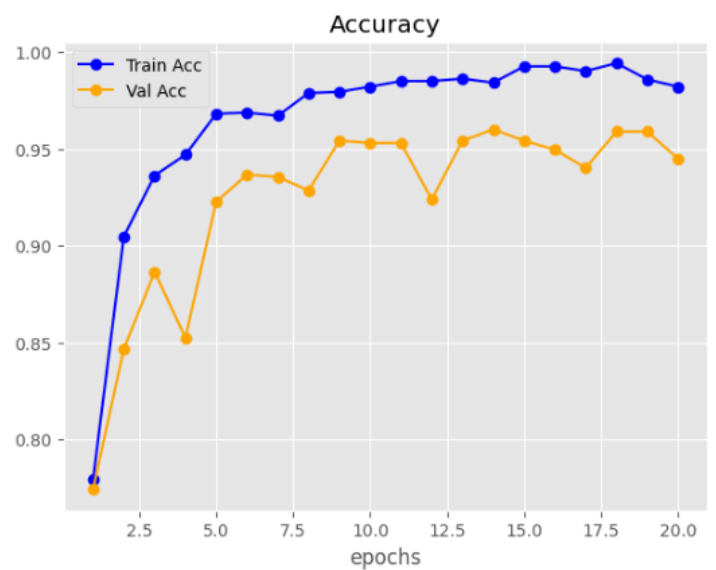


FIGURE 9: ACCURACY GRAPH (VGG16)

### 3.6. RESNET101 MODEL ARCHITECTURE:

To further enhance the classification performance and improve model generalization, a deeper convolutional neural network architecture, ResNet101, was employed in addition to the VGG16 model. Residual Networks (ResNet) are specifically designed to address the limitations of very deep neural networks, particularly the degradation and vanishing gradient problems that occur as network depth increases.

ResNet introduces the concept of **residual learning**, where shortcut (skip) connections allow the input of a layer to bypass one or more intermediate layers and be added directly to the output. Instead of learning a direct mapping, the network learns a residual function. This design enables the model to train effectively even with a large number of layers.

Mathematically, the residual block can be expressed as:

$$H(x) = F(x) + x$$

where:

- $x$  represents the input,
- $F(x)$  is the learned residual mapping,
- $H(x)$  is the final output of the block.

These shortcut connections help preserve gradient flow during backpropagation, allowing deeper networks to converge faster and achieve higher accuracy.

ResNet101 consists of **101 layers**, including multiple convolutional layers organized into residual blocks, batch normalization layers, and activation functions. Compared to VGG16, ResNet101 is significantly deeper and capable of extracting more complex and discriminative features from MRI images.

In this project, ResNet101 was trained using the merged dataset that combines the original MRI images with newly collected samples. The deeper architecture allows the model to learn richer representations and improve robustness when handling variations in tumor shapes, sizes, and image conditions.

The final fully connected layer of the network was modified to match the number of tumor classes in the dataset, enabling the model to perform multi-class brain tumor classification.

## Base Model Configuration

- Pretrained Weights: ImageNet
  - Input Shape: 150×150×3 (RGB MRI slices)
  - Frozen Layers: All layers except the last 20 layers
  - Trainable Layers: Last 20 layers for fine-tuning
  - Regularization: 50% Dropout
- 

## Key Modifications from Original ResNet101

- Removed original 1000-class fully connected layer (`include_top=False`)
  - Added Global Average Pooling layer
  - Added Dense layer with 512 neurons (ReLU activation)
  - Added Dropout layer for overfitting reduction
  - Replaced final classifier with Dense layer of *num\_classes* outputs (Softmax)
- 

## Training

- Optimizer: Adam (`learning_rate = 0.0001`)
  - Loss Function: Categorical cross-entropy
  - Metrics: Accuracy
-

## Training Protocol

- Epochs: 20
- Batch Size: (same as data generator setting 32)
- Validation Data: validation generator
- Transfer Learning Strategy: freeze early layers and fine-tune final layers only

## Confusion Matrix

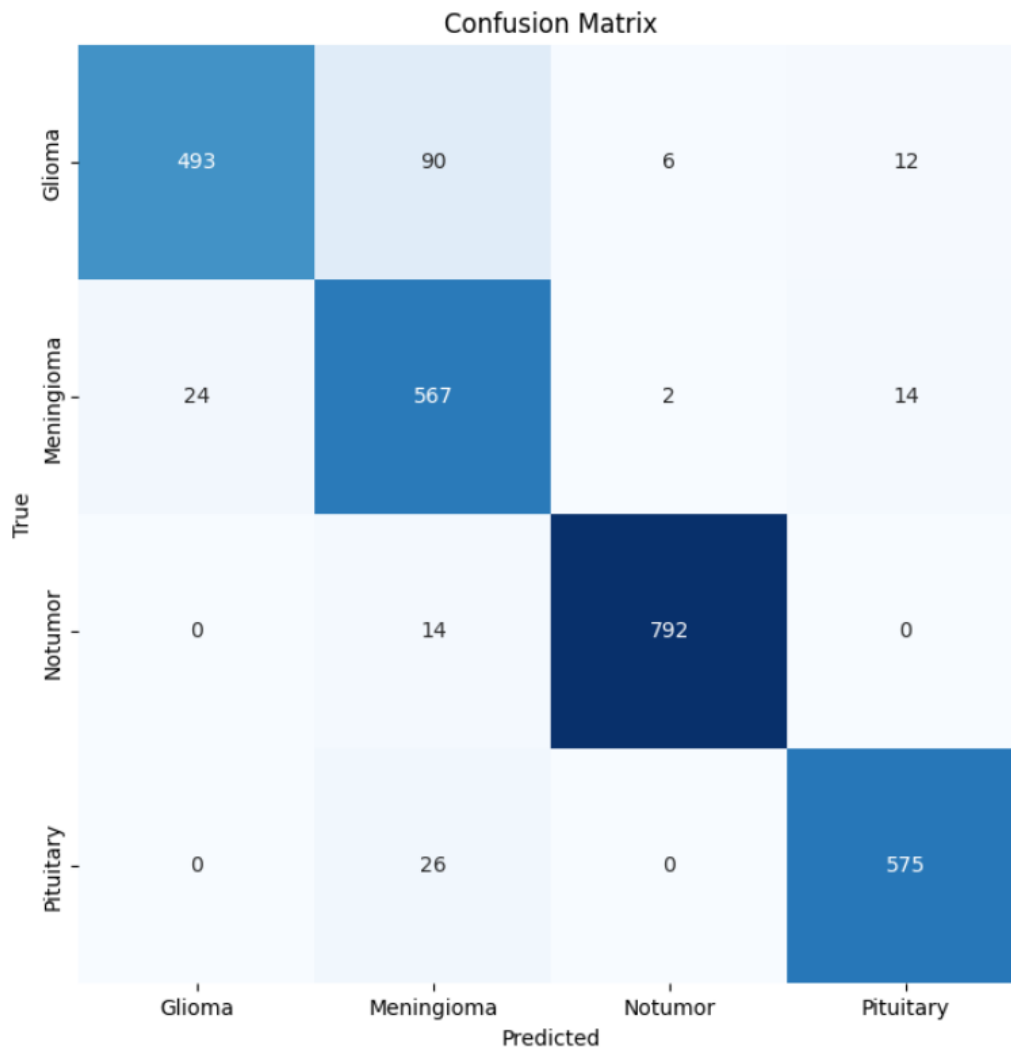


FIGURE 10: CONFUSION MATRIX (RESNET101)

## Accuracy & Loss (RESNET101)

```
loss, accuracy = model_resnet.evaluate(test_generator)
print(f"Test Loss: {loss:0.5f}")
print(f"Test Accuracy: {accuracy:0.5f}")
```

**82/82** ————— **58s** 704ms/step - accuracy: 0.8773 - loss: 0.3644  
Test Loss: 0.21205  
Test Accuracy: 0.92811

FIGURE 11 TEST ACCURACY & LOSS (RESNET101)



FIGURE 13 LOSS GRAPH (RESNET101)

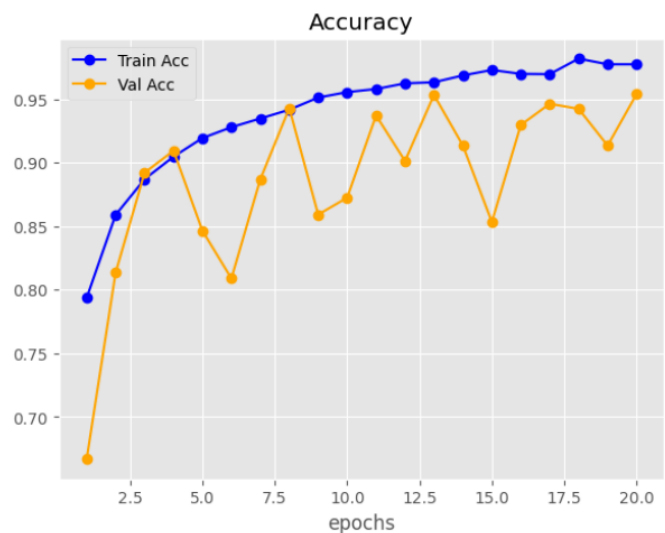


FIGURE 12 ACCURACY GRAPH (RESNET101)

## 4. Chapter 4: Android application implementation

### 3.7. System analysis & Design:

#### 3.7.1. Requirements

TABLE 4-1 REQUIREMENTS

ID	Requirement Type	Description	Priority	Remarks
FR-01	Functional	User can upload MRI images from device storage	High	Single image at a time
FR-02	Functional	System preprocesses uploaded image (resize 150x150 + normalize)	High	Ensures consistent input for models
FR-03	Functional	System runs VGG16 and ResNet101 models on image	High	TFLite float16 models
FR-04	Functional	System compares model predictions; if both agree, display tumor type and average confidence	Medium	Voting ensemble mechanism
FR-05	Functional	If model predictions conflict, show “unreliable scan” message	Medium	Prompt user to provide clearer MRI
FR-06	Functional	Display individual model outputs in debug box	Medium	Helps in debugging or confidence explanation
FR-07	Functional	Vibrate device if tumor is detected	Low	Optional UX enhancement
FR-08	Functional	Provide progress bar during inference	Low	Shows processing status
FR-09	Functional	Help button toggles debug info visibility	Medium	Improves user guidance
NFR-01	Non-Functional	Models must be loaded asynchronously	High	Prevent UI blocking
NFR-02	Non-Functional	Support devices compatible with Android 11+	High	Minimum Android API 30 recommended
NFR-03	Non-Functional	Predictions must complete within 3 seconds for standard images	Medium	Dependent on device hardware
NFR-04	Non-Functional	UI must be responsive and follow Material Design	High	Uses ConstraintLayout and Material components
NFR-05	Non-Functional	TFLite models must not be compressed in APK	High	Ensures models load correctly



### 3.7.2. Use Case Diagram:

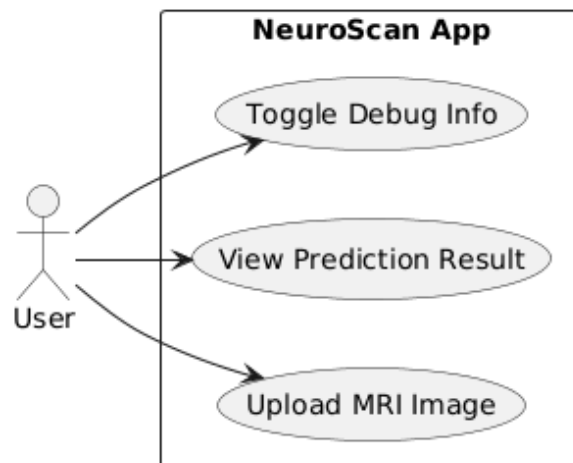


FIGURE 14 USE CASE DIAGRAM

### 3.7.3. Activity Diagram

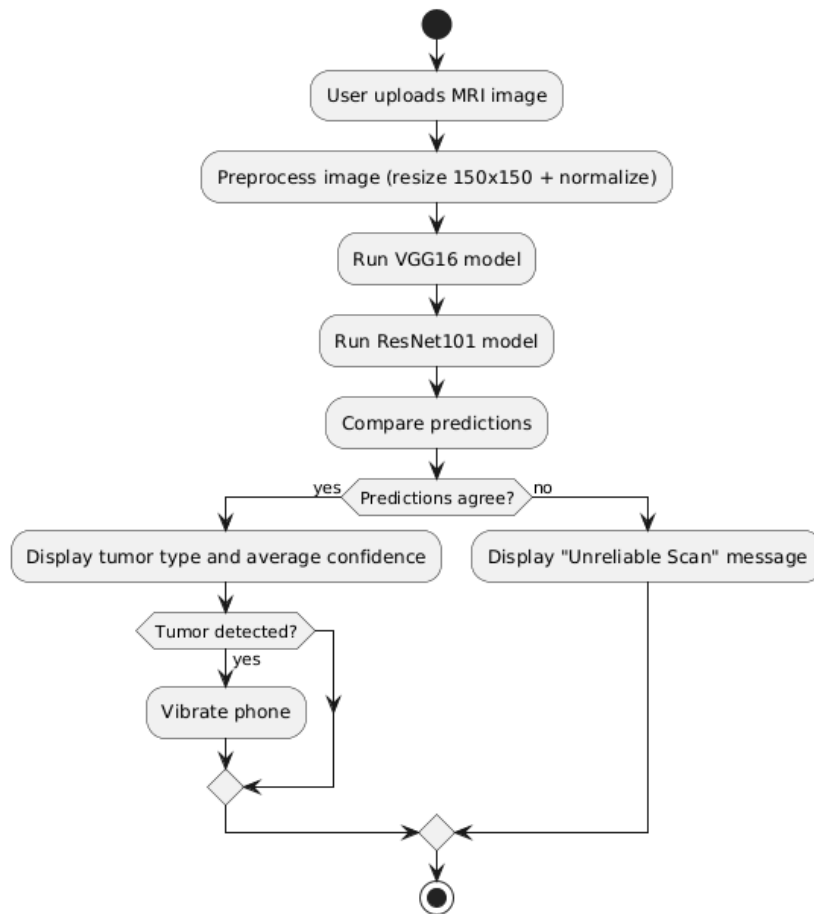


FIGURE 15 ACTIVITY DIAGRAM

### 3.7.4. Sequence Diagram

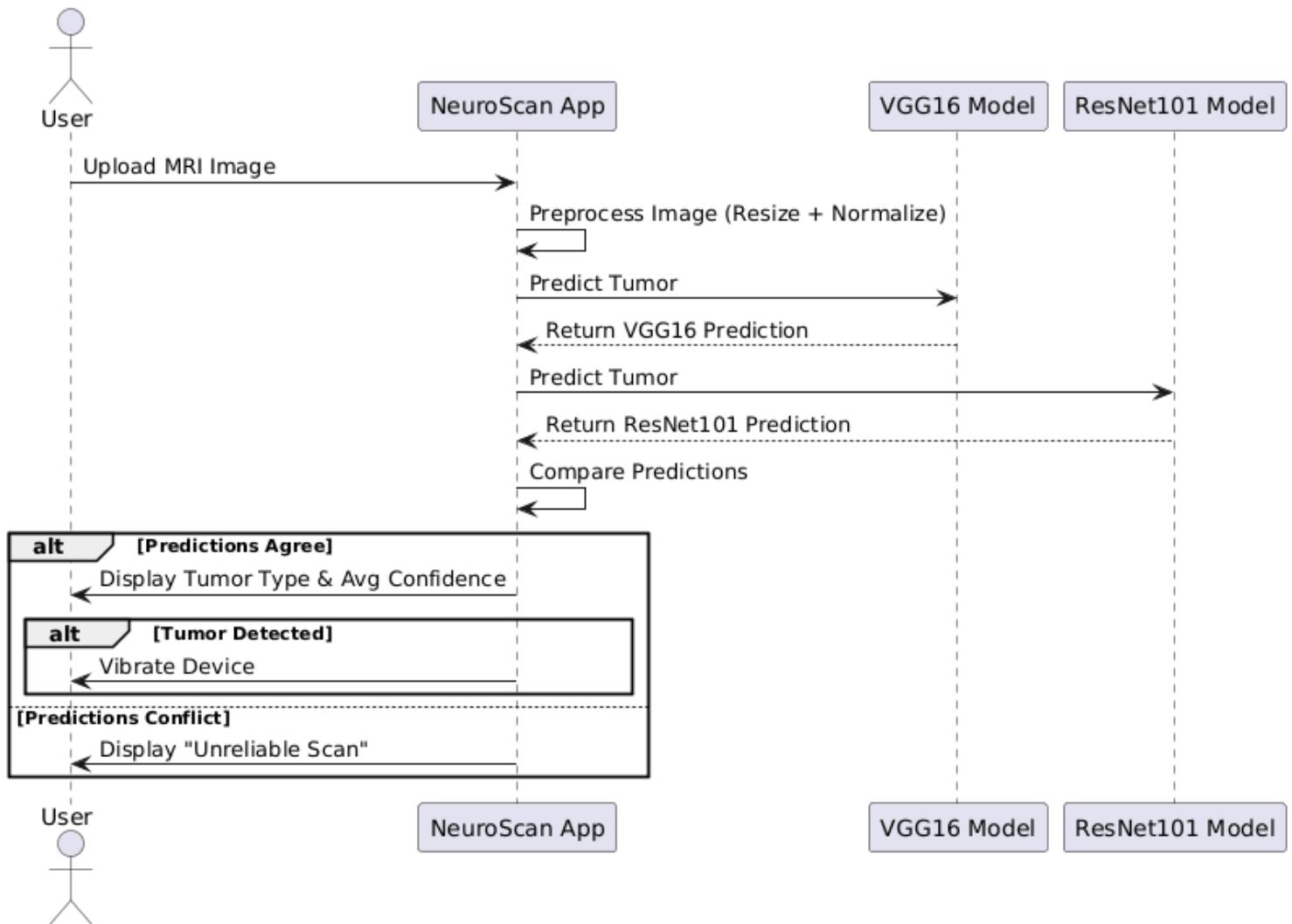


FIGURE 16 SEQUENCE DIAGRAM

### 3.7.5. General System Schema (High-Level Architecture)

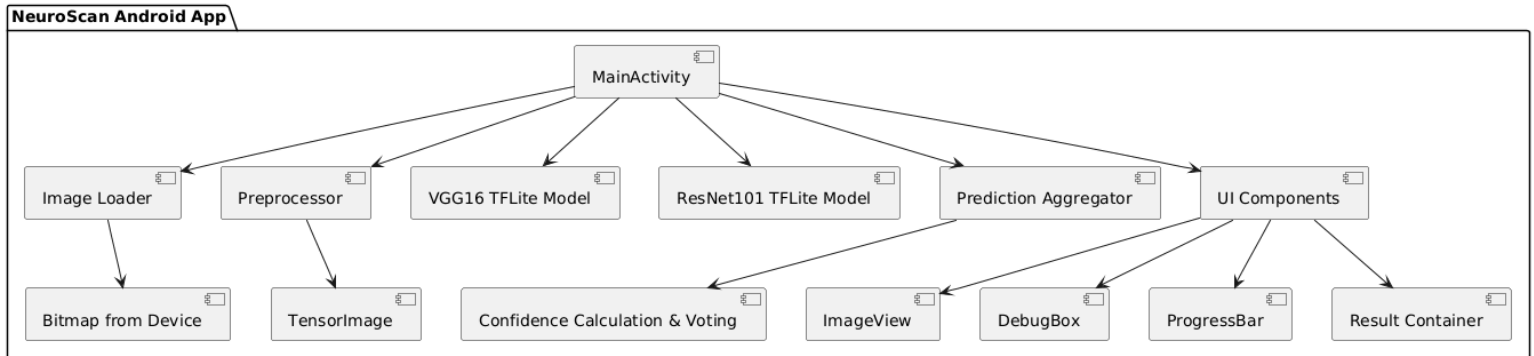


FIGURE 17 GENERAL SYSTEM SCHEMA

### 3.7.6. Class Diagram

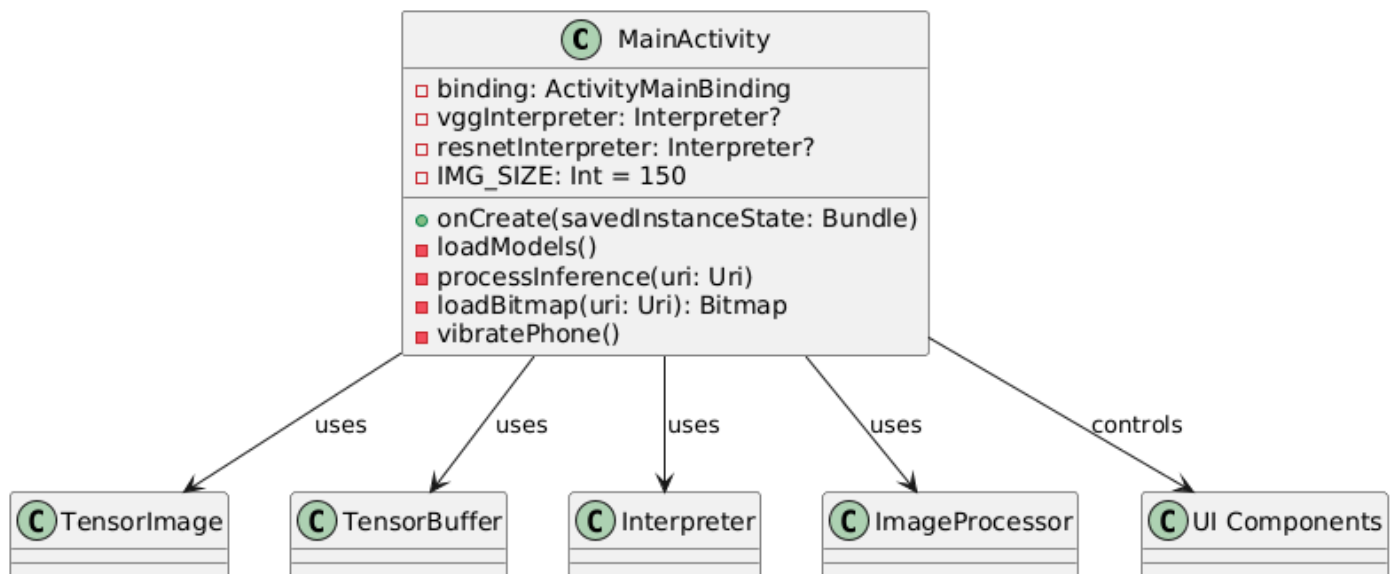


FIGURE 18 CLASS DIAGRAM

#### 4.1.7. Test Cases

**TABLE 4-2 TEST CASES**

Test Case ID	Test Scenario	Test Steps	Input Data	Expected Result	Priority
TC-01	Upload MRI Image	1. Launch app 2. Tap “Analyze New Scan” 3. Select MRI image	Valid MRI image	Image displayed in ImageView	High
TC-02	Preprocess Image	1. Upload MRI image 2. Preprocessing triggered automatically	Valid MRI image	Image resized to 150x150, normalized	High
TC-03	VGG16 Prediction	1. Preprocess image 2. Run VGG16 inference	Valid MRI image	Model outputs 4-class probabilities	High
TC-04	ResNet101 Prediction	1. Preprocess image 2. Run ResNet101 inference	Valid MRI image	Model outputs 4-class probabilities	High
TC-05	Ensemble Voting — Agree	1. Run both models on image 2. Both models predict same class	Valid MRI image	Display tumor type & average confidence	High
TC-06	Ensemble Voting — Conflict	1. Run both models on image 2. Models predict different classes	Valid MRI image	Display “Unreliable Scan” message	High
TC-07	Confidence Display	1. Models agree 2. Check result box	Valid MRI image	Confidence shown in percentage (avg of both models)	High
TC-08	Vibrate on Tumor	1. Tumor detected 2. Check device response	MRI with tumor	Device vibrates 300ms	Medium
TC-09	No Vibrate for No Tumor	1. “No Tumor” detected	MRI with no tumor	Device does NOT vibrate	Medium
TC-10	Debug Box Toggle	1. Tap Help FAB 2. Tap again	Any	Debug box visibility toggles	Medium
TC-11	Progress Bar	1. Upload image 2. During inference	Any MRI	Progress bar visible while processing, hides when done	High

Test Case ID	Test Scenario	Test Steps	Input Data	Expected Result	Priority
TC-12	Invalid File Type	1. Try to upload non-image file	PDF / TXT / etc.	App shows error or fails gracefully	High
TC-13	Large Image	1. Upload high-res MRI (>5MB)	Large MRI image	Image processed correctly without crash	High
TC-14	Multiple Sequential Uploads	1. Upload image 2. After result, upload another image	Multiple MRI images	Each image processed independently	Medium
TC-15	System Crash/Exception	1. Force failure in model loading	Corrupt TFLite file	App shows error message, does not crash	High

## Chapter 5: Conclusion & Recommendations

## 3.8. Conclusion

In this chapter, the implementation of the Android application was presented in detail. The chapter covered the system architecture, main modules, user interface design, and integration of dual deep learning models (VGG16 and ResNet101) through a voting-based ensemble mechanism. The workflows, from image upload and preprocessing to inference, confidence calculation, and result presentation, were illustrated using UML diagrams including use case, activity, and sequence diagrams.

The chapter also addressed the app's robustness through handling prediction conflicts, error messages, and edge cases such as invalid or large MRI images. Functional and non-functional requirements were mapped to specific features and test cases to ensure reliability and usability. Furthermore, performance considerations such as asynchronous model loading and use of float16 models were incorporated to maintain a smooth and responsive user experience.

Overall, the implementation demonstrates a seamless integration of AI models into a mobile application while maintaining a user-friendly interface and adherence to software engineering best practices. The NeuroScan app successfully provides an efficient, reliable, and interpretable platform for brain tumor detection on Android devices.

### 3.9. Recommendations

Based on the outcomes of this project, several recommendations are proposed to further enhance the system and extend its capabilities:

- Increasing the dataset size with more diverse MRI scans to improve model robustness and generalization
- Evaluating additional deep learning architectures to compare accuracy and computational efficiency
- Applying advanced data augmentation techniques to reduce overfitting
- Optimizing model size and inference speed for lower-end mobile devices
- Integrating direct camera capture in addition to gallery image selection
- Adding a history or logging feature to store previous predictions for review
- Providing detailed probability distributions for all tumor classes
- Improving the user interface for better usability and accessibility
- Conducting clinical validation in collaboration with medical professionals to assess real-world reliability
- Extending the system to cloud or web-based platforms for remote access and large-scale deployment

These improvements can enhance system performance, usability, and practical applicability in medical and research environments.



## References:

1. Abdusalomov, Akmalbek Bobomirzaevich, Mukhriddin Mukhiddinov, and Taeg Keun Whangbo. "Brain tumor detection based on deep learning approaches and magnetic resonance imaging." *Cancers* 15.16 (2023): 4172.
2. Srinivas, Chetana, et al. "Deep transfer learning approaches in performance analysis of brain tumor classification using MRI images." *Journal of Healthcare Engineering* 2022.1 (2022): 3264367.
3. Mokri, S. M., Newsha Valadbeygi, and Vera Grigoryeva. "Diagnosis of Glioma, Meningioma and Pituitary brain tumor using MRI images recognition by Deep learning in Python." *EAI Endorsed Transactions on Intelligent Systems and Machine Learning Applications* (2021).
4. Haq, Amin ul, et al. "DACBT: Deep learning approach for classification of brain tumors using MRI data in IoT healthcare environment." *Scientific Reports* 12.1 (2022): 15331.
5. Nalepa, Jakub, et al. "Fully-automated deep learning-powered system for DCE-MRI analysis of brain tumors." *Artificial intelligence in medicine* 102 (2020): 101769.
6. Anaraki, Amin Kabir, Moosa Ayati, and Foad Kazemi. "Magnetic resonance imaging-based brain tumor grades classification and grading via convolutional neural networks and genetic algorithms." *biocybernetics and biomedical engineering* 39.1 (2019): 63-74.
7. Yousef Mohamed, 2024. Brain Tumor MRI Datase, A dataset for classify brain tumors. Brain tumor MRI

**Note:**

.....

.....

.....

.....

.....

.....

.....

.....

.....